


826

Full Face Mask Nasal Mask

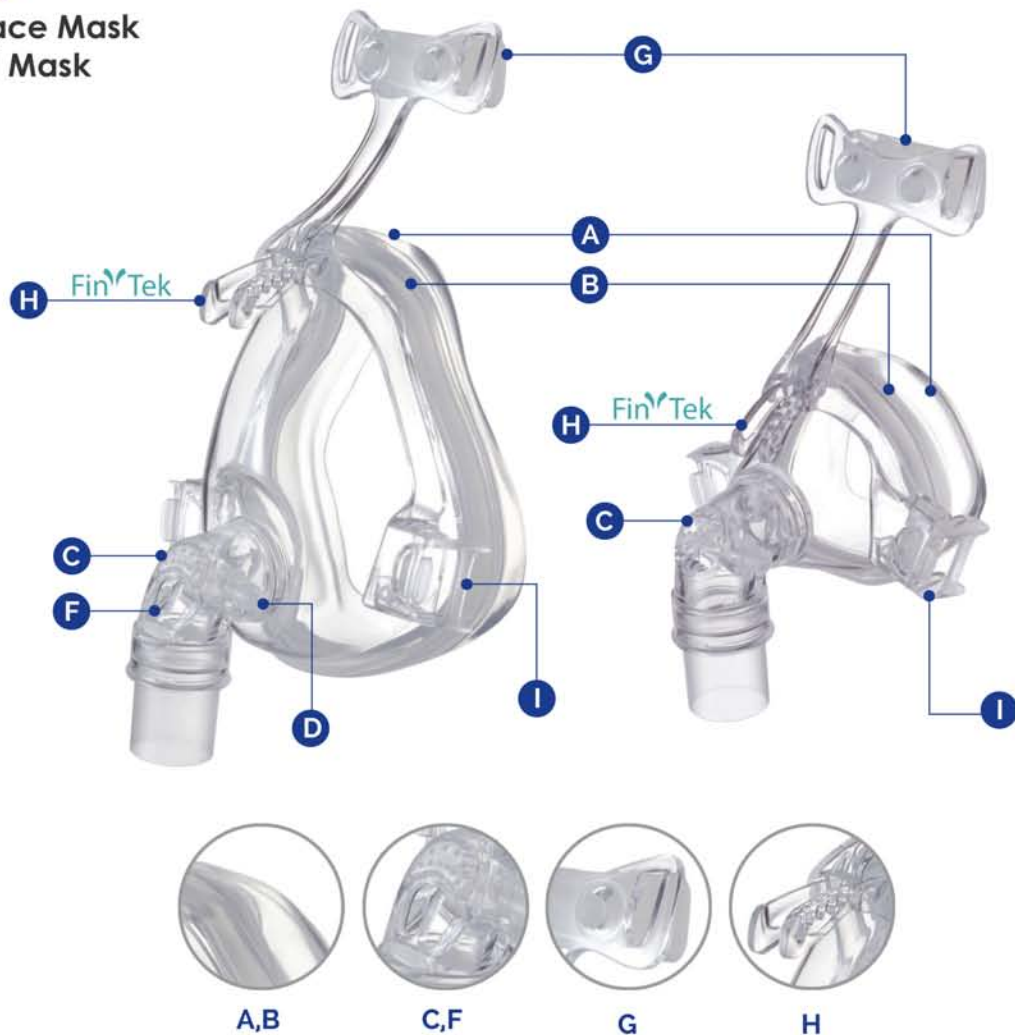

- A** Medical grade liquid silicone full face cushion fits just right and firmly on your face.
- B** Quick assembling design offer user more convenience when taking off / assembling the cushion.
- C** Silence type vent holes provides lower noise and the maximum CO₂ washout, vented air is spread as cloud so that does not disturb patient or companion.
- D** Sampling port allows measuring pressure or ETCO₂.
- E** Minimum components, less components will lower the risk of losing components and easier for user when cleaning.
- F** Anti-asphyxia valve, automatically open when the generator pressure fails, enabling patients to breathe by themselves. (Full Face Mask)
- G** Flexible, One-piece forehead pad evenly conforms to the forehead contour, so the pad stays in place and you won't feel a thing.
- H** **FinTek** >>Innovative adjustable forehead support provides stability with easy one-handed adjustment and allows the mask fit most of patients.
- I** Push release / attach buckle, easier for use to take off mask after sleep.