


837

Pediatric Total Face Mask

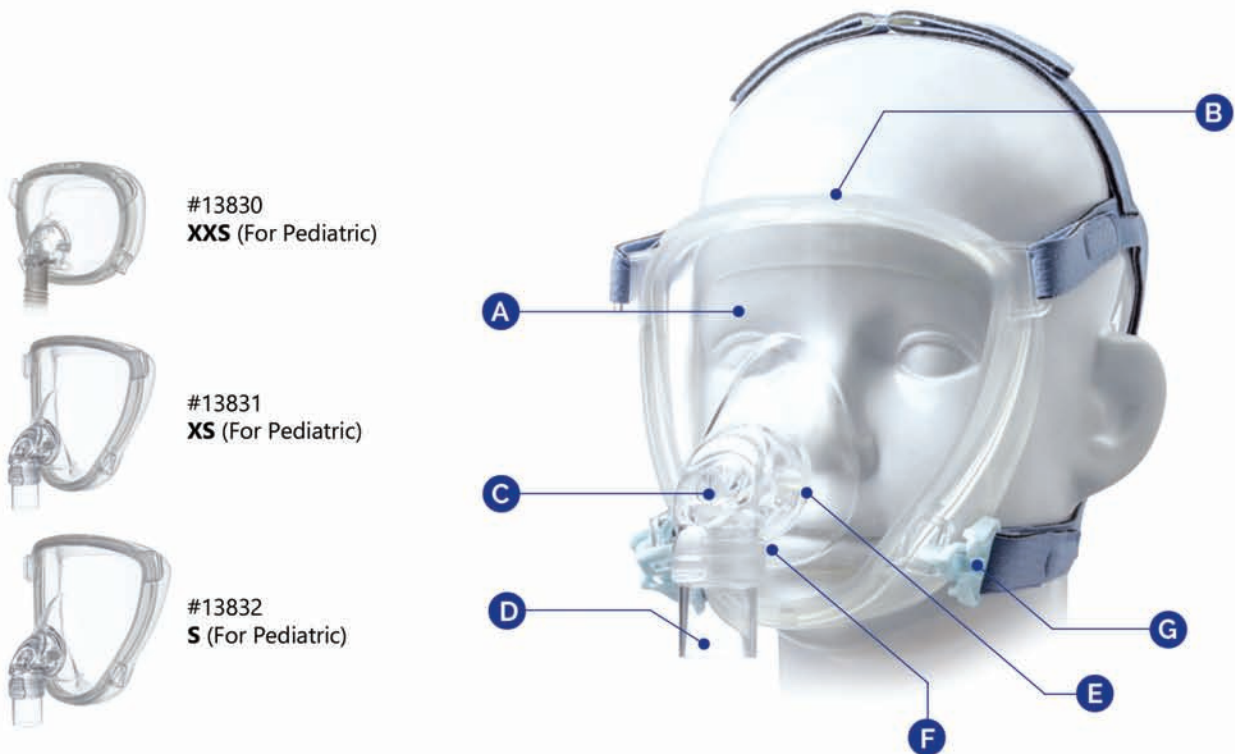

- A** The transparent cover creates a crystal clear view for the user.
- B** Silicone total full face cushion fits just right and firmly on the whole face, has less pressure for the user.
- C** Anti-asphyxia valve, automatically open when the generator pressure fails, enabling patients to breathe by themselves
- D** Quick release adaptor, press make user release / attach the tubing easily during temporary leaving.
- E** Sampling port allows measuring pressure or ETCO₂.
- F** With nasogastric (NG) port.
- G** Quick release buckles, easy for user to release and take the mask after sleep.
- H** 360° rotational elbow providing a comfortable tube position.
- I** Applicable to the users who have facial hair and general facial irregularity.